



Rumble Tuff®

Rumble Tuff
Electric Breast Pump
Instructions & Assembly
for Mighty Tyke

Congratulations

Congratulations on your new Rumble Tuff® Electric Breast Pump! The Rumble Tuff® Electric Breast Pump is perfect for many occasions.

- Collecting breast milk while you're away from your baby
 - for example: when you're at work
- Relieving engorgement
 - for example: when your baby's feeding habits change because they're eating more solids
- Maintaining your milk supply if you have to temporarily stop breast feeding
 - for example: to express unsuitable milk while you are undergoing medical treatment

The Rumble Tuff® Electric Breast Pump provides many of the features busy moms love.

- 2 cycles of expressions to mimic baby's nursing patterns.
 - Swift Expression - Fast but shallow swallowing habit
 - Natural Expression - Steady and strong swallowing pattern
- Milk back-flow protection
- Quiet motor for discretion
- Small and lightweight
- 2 power options to fit your needs
 - AC/DC Adapter
 - Li-ion battery

Indication for Use

The Rumble Tuff® Electric Breast Pump is a single user, electrical powered device used to express milk from the breasts of lactating women. This device is not intended for hospital use.

Table of Contents

Congratulations	0
Indication for Use	0
Important Safeguards	3
• DANGER	3
• CAUTION	3
• WARNING	3
Parts Guide	5
• Electrical Pump Main Unit	6
• LED Screen Guide	6
• Expression Collection Combo Kit	7
How to Assemble the Expression Collection Kit	8
How to Assemble Electric-pump Diaphragm Kit	11
How to Assemble Expression Collection Combo Kit	12
Power Options	14
• Using your AC/DC Adapter	14
• Using the Li-Ion Rechargeable Battery	15
• Install the Battery	16
Pumping Instructions	18
LED Light Signal Guide	22
Cleaning Instructions	23
• Sanitize Prior to First Use	23
• Before Each Use	23
• General Cleaning Instructions	23
• Tubing Cleaning Instructions	24
Breastmilk Storage Guidelines	25
• Important Notes	25

Tips	27
• Pumping Tips	27
• Battery Storage	28
• Breastfeeding Tips	28
Troubleshooting	29
• Interruption of Electrical Service	29
• Radio Frequency Energy	29
• Decreased Suction	29
• Pumping Doesn't Start	30
• Breastmilk Buildup in the Breastshield Body	30
• Where to get Replacement Parts	30
• Where to get Additional Accessories	30
• What Additional Accessories are Available	30
Technical Data	31
• Power Requirements	31
Warranty	32

Important Safeguards

The Rumble Tuff® Electric Breast Pump is designed for a single user and should never be shared or resold; doing so will void the warranty.

READ ALL INSTRUCTIONS BEFORE USING

DANGER: To reduce the risk of electrocution:

- Always unplug this product immediately after using.
- Do not use while bathing.
- Do not place or store this product where it can fall or be pulled into a tub or sink.
- Do not place or drop into water or other liquid.
- Do not reach for an electrical product that has fallen into water. Unplug immediately.

CAUTION:

- Do not use this product for any other purpose than as indicated in the manual.
- Make sure all the parts are assembled correctly before operating.
- Stop using the breast pump immediately if your breasts or nipples experience any discomfort.
- Please wash hands before using this product.
- Keep this product out of reach of babies and children.
- Do not drop or hit the product to avoid damage.
- Do not touch the product with bare hands right after sterilizing by boiling water to avoid scalds and burns.
- Do not over-boil to avoid premature aging of the product materials.

WARNING: To reduce the risk of burns, electrocution, fire, or injury to persons:

- This product should never be left unattended when plugged in.
- Close supervision is necessary when this product is used near children or invalids.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by Rumble Tuff®.
- Do not operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it has been dropped in water.
- Keep the cord away from heated surfaces.
- Do not use while sleeping or drowsy.
- Do not drop or insert any other object into any opening or tubing.
- Do not use outdoors with the AC/DC Adapter or operate where aerosol spray products are being used, or where oxygen is being administered.

SAVE THESE INSTRUCTIONS
along with your receipt for future reference.


Important Safeguards (Continued)

- Make sure the voltage of the AC/DC Adapter is compatible with the power source.
- Plug the AC/DC Adapter into the electrical outlet first, then plug the Power Cord into your Rumble Tuff® Electric Breast Pump.
- Rumble Tuff has provided an AC/DC Adapter and Li-Ion battery pack with this device. Using any other AC/DC Adapter or Li-Ion battery pack not approved by the manufacturer voids the warranty.
- If you are infected with Hepatitis A, Hepatitis B, or have HIV (Human Immunodeficiency Virus): pumping breast milk does not reduce or remove the risk of transmitting the virus to your baby through breast milk.
- Wash the Expression Collection Combo Kit parts before each use.
- Do not operate a motor vehicle while using the Rumble Tuff® Electric Breast Pump.
- Do not attempt to modify the product.

Parts Guide




Mighty Tyke Breast Pump Main Unit

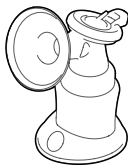


Expression Collection Combo Kit



Vacuum Tube with Adapter

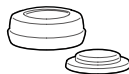


Bottle Adapter Kit

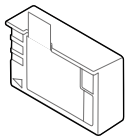

Li-Ion
Rechargeable Battery




Extra 5 oz. Collection Bottle

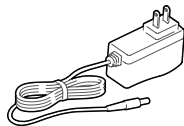


AC/DC Adapter

Parts Guide

• Electric Pump Main Unit

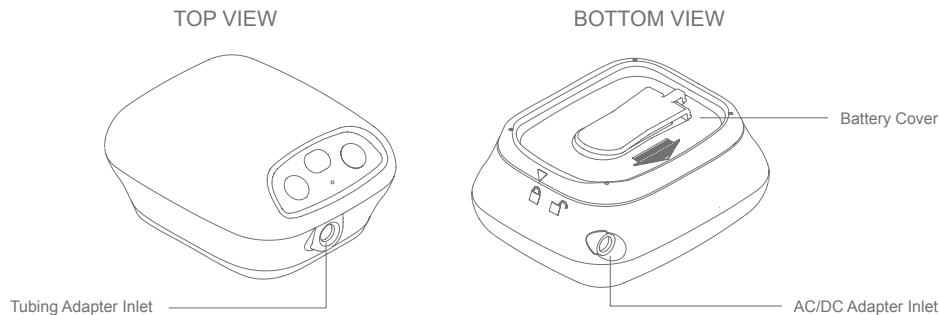


• LED Screen Guide

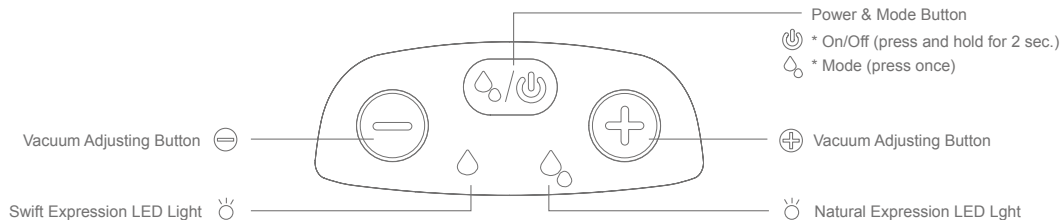

Parts Guide

● Expression Collection Combo Kit

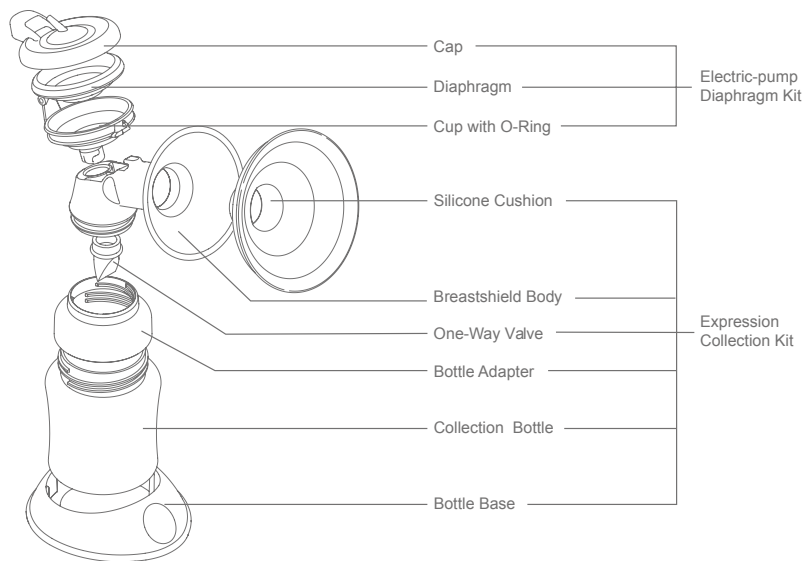

How to Assemble the Expression Collection Kit


Your Rumble Tuff® Electric Breast Pump comes pre-assembled; however, you will need to re-assemble it after each cleaning. (For cleaning instructions, see page 23) Below is a step-by-step guide for assembling the Expression Collection Kit:

1. Hold the "horn" of the Breastshield Body in one hand. Make sure the bottle-thread side opening is facing up.



2. With the thumb and forefinger of the other hand, hold the tip of the One-Way Valve and place the flat circular bottom into the opening. Be sure not to stretch, bend, or twist the tip of the One-Way Valve. Seal all the edges firmly.


Note: The One-Way Valve tip does not have to align with anything, as long as all the edges are securely sealed.

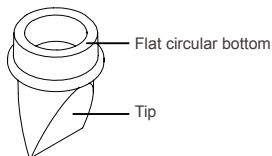

How to Assemble the Expression Collection Kit (Continued)

3. With your thumb and forefinger, gently seal any openings on the tip of the One-Way Valve and shape into a straight line.

Note: This is critical to the performance of the Pump.

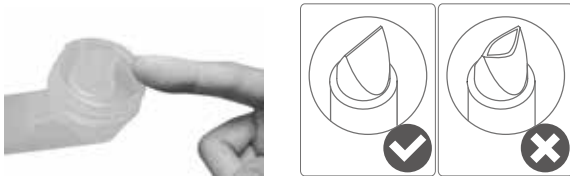

4. Carefully twist on the Bottle Adapter to the Breastshield Body for a secure fit.

Note: Be sure not to bump the One-Way Valve when twisting the Breastshield Body & the Bottle Adapter pieces together. If you bump the One-Way Valve, repeat step 3.



How to Assemble the Expression Collection Kit (Continued)

5. With the Bottle Adapter in one hand, twist on the Collection Bottle for a secure fit.



How to Assemble Electric-pump Diaphragm Kit

1. With the Cup in one hand, place the Silicone Diaphragm in the Cup (matching the shape of the Cup). Press around the edges of the Diaphragm for a secure fit. Make sure the Diaphragm is evenly mounted on the Cup and there is no gap in between them.




2. Place the Cap on the Cup & Diaphragm. (Be sure to align the hinge and hinge cover.) Press tightly to seal.



How to Assemble Expression Collection Combo Kit

1. With the O-Ring side down, put the Electric Pump Diaphragm Kit into the top of the Breastshield Body opening. Twist the Electric-Pump Diaphragm Kit clockwise into the Breastshield Body until it locks into place.

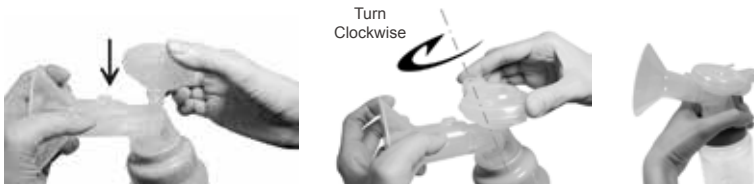

2. With the Breastshield Body in one hand, place the Silicone Cushion in the Breastshield Body (matching the shape of the horn). Secure the lip of the Cushion around the edges of the horn.


Note: Your Electric Breast Pump will work with or without the Silicone Cushion. Your comfort depends on your preference.



How to Assemble Expression Collection Combo Kit (Continued)

3. Take the open-end of Vacuum Tube (the end without the adapter) and connect to the Tubing Connector located on the Cap. Push together firmly for a tight seal.

Note: Make sure there is no moisture in the tubing.

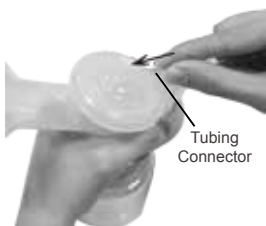

4. Make sure the Tubing Adapter is properly inserted into the Tubing Adapter Inlet on the Electric Pump Main Unit.



Power Options

Your Electric Breast Pump may be operated under any of the following power supply modes:

Using your AC/DC Adapter:

- Electrical safety can only be achieved when using the AC/DC Adapter which was provided with the Pump. Use of any another Adapters causes concern for the safety of the product and/or the user and voids the warranty.


1. Plug the AC/DC Adapter into an electrical outlet.



2. Attach the Power Cord to the AC/DC Adapter Inlet on the Electric Pump Main Unit.



Using the Li-Ion Rechargeable Battery:

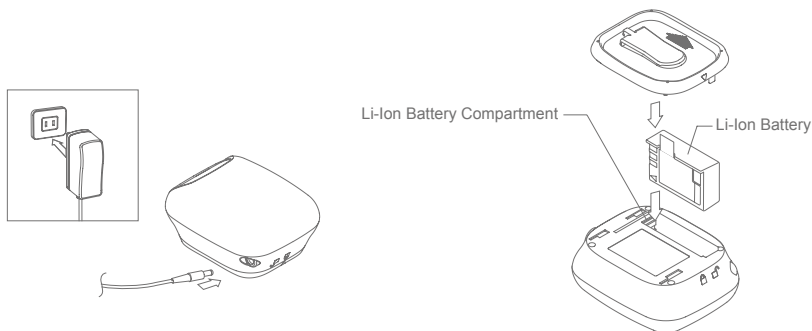

- Prior to first use, charge the battery for 24 hours. If necessary, your Electric Breast Pump can be used during this period with the Power Cord only. For best results, allow the Pump to initially charge uninterrupted for 24 hours.
- After the initial charge, subsequent charges will take approximately 4 hours. Once charged, you should not need to remove the battery unless replacing or when storing.
- The Li-Ion Rechargeable Battery will charge every time the AC/DC Adapter is plugged in; however, it will not charge while your Electric Breast Pump is in use.
- The Swift Expression LED and Natural Expression LED flash alternately while the battery is charging.

Note: *The Li-Ion Rechargeable Battery used in this device may present a risk of fire or chemical burn if mistreated. Do NOT disassemble the battery. Do NOT place battery near heat source or fire. Dispose of used battery properly. Do NOT dispose of battery in fire. For long term storage, keep battery within the following temperature range: 32°F to 77°F (0°C to 25°C). Keep away from children.*


Install the Battery: Li-Ion Rechargeable

1. Open the Battery compartment by sliding the cover to the side by following the arrow. Lift to remove the cover.



2. Use the tab to remove the old battery.

Note: To remove the battery, slide the battery down, away from the side with the Tubing Adapter Inlet. Pull the tab up to remove the old battery.

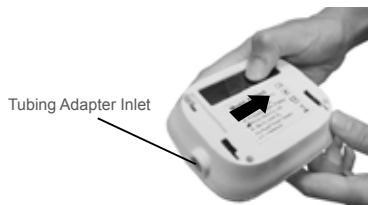

Install the Battery: Li-Ion Rechargeable (Continued)

3. Insert the new Li-Ion Rechargeable Battery. Make sure the sticker on the battery faces the sticker on the pump's bottom, and the battery is inserted completely into the compartment.

Note: The pump will not operate on battery power if the battery is inserted incorrectly. When the battery is inserted correctly, it will slide easily into the battery compartment. Do NOT force.




4. Close the Battery Cover and slide the back in place until locked. Align the arrow to the locked position to ensure the battery cover is closed.

Note: Make sure the battery is inserted completely into the compartment before pushing the cover closed.



Pumping Instructions


1. Hold the Expression Collection Combo Kit to your breast so that your nipple is centered inside the Breastshield.

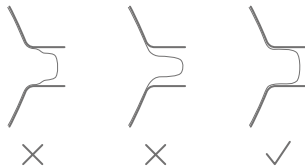

Note: Please check for correct fitting of the shield.

2. Press and hold the Power/Mode Button  until the pumping starts.


Note: Make sure to hold the Power/Mode button long enough (approx. 2 sec.) until the pumping starts. If this doesn't work, please check the Troubleshooting section on page 29.

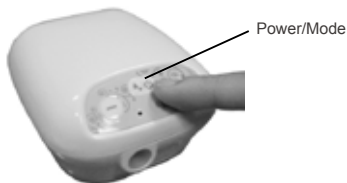



Pumping Instructions (Continued)

3. The Pump will start in the Swift Expression mode  at the lowest suction level. The Swift Expression LED will light up.

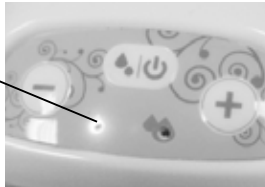Swift Expression LED Light



4. Adjust the suction to your comfort level while pumping by pressing the buttons  or . A beep sound can be heard whenever the suction level is changed.

Note: The LED light will flash when you have reached the minimum or maximum vacuum settings (There are 8 suction levels).

Vacuum Adjusting buttons



Pumping Instructions (Continued)

5. If you feel your milk start to let-down, press the Power/Mode Button  to switch to the Natural Expression mode. 

Note: The pump will turn off if you press and hold the Power/Mode Button for more than 2 seconds.

Power/Mode Button



6. The pump will switch to the Natural Expression mode at the lowest suction level. The Natural Expression LED will light up.

Swift Expression mode



Natural Expression mode



Pumping Instructions (Continued)

7. Adjust the suction in the Natural Expression mode to your comfort level while pumping by pressing the Vacuum Adjusting buttons



Note: The LED light will flash when you have reached the minimum or maximum vacuum settings (There are 8 suction levels).

Vacuum Adjusting
buttons



8. When your pumping session is finished, press and hold the Power/Mode Button  for more than 2 seconds until the pumping stops.

Note: Make sure to hold the Power/Mode button long enough (approx. 2sec.) until the pumping stops.
















Make sure there is no suction on the nipple before removing the Breastshield. If necessary, break the suction by inserting a finger between the Breastshield and the breast, similar to the way you break away from your baby's suction.

Power/Mode



Note: If the Rumble Tuff® Electric Breast Pump is left on, but unattended for more than 30 minutes, the Pump will automatically shut-off.

LED Light Sign Guide

	State	Signal
While Pump is on	Pumping in Swift Expression	<div>Swift Expression LED (White light)</div>  <div>White Light is On</div> <div>Natural Expression LED (Blue light)</div> 
	Pumping in Natural Expression	<div>Blue Light is On</div>  
	Battery is Low (Swift Expression mode)	<div>White Light Flashes</div>  
	Battery is Low (Natural Expression mode)	<div>Blue Light Flashes</div>  
	Battery is Empty	<div>Both Blue & White Lights Flash, then turn Off</div>  
While Pump is off	Battery is Charging	<div>Blue & White Lights Flash Alternately</div>   
	Charging is Completed	<div>Both Blue & White Lights are On</div>  

 Light is Off
  Light is On
  Light is Flashing

Cleaning Instructions

Unless otherwise instructed by your healthcare provider, follow these sanitizing and cleaning instructions.

Sanitize Prior to First Use:

Before you use your Rumble Tuff® Electric Breast Pump for the first time, disassemble and boil all parts of the Expression Collection Combo Kit in water for 10 minutes.

- Be sure to use enough water to keep all parts afloat while boiling.
- **Do NOT boil Tubing or Tubing Adapter.**
- Remove all parts from the water immediately after sterilizing. To avoid scalds and burns, do NOT touch the product with bare hands right after sterilizing.
- Completely air-dry or towel-dry before reassembling and next use.

Partially dry components may affect the operation of your Electric Breast Pump. Please ensure that all parts are completely dry before use.

Before Each Use:

- Wash your hands thoroughly with soap and warm water before handling the Expression Collection Combo Kit, the Electric Breast Pump or any Expression Collection parts or attachments.
- Make sure all parts are clean and dry.
- Ensure the Tubing is free from any liquid. **Do NOT operate your Rumble Tuff Electric Breast Pump if any liquid is visible in the Tubing. (See Tubing Cleaning Instructions on page 24.)**

General Cleaning Instructions:



Caution: Take care when cleaning the One-Way Valve. If it gets damaged, your breast pump will not function properly. To clean the One-Way Valve, rub it gently between your fingers in warm soapy water, then rinse thoroughly in clean warm water. Do not insert any objects into the One-Way Valve as this may damage the valve.

- Set aside the Tubing and Tubing Adapter. Under the normal working conditions of the pump, the diaphragm prevents the back-flow of milk. Therefore, these items do not need to be washed. In case of finding any liquid or liquid vapor inside the Tubing, follow the Tubing Cleaning Instructions (page 24) carefully before use.
- All parts of the Expression Collection Combo Kit may be cleaned in the upper compartment of your dishwasher or hand washed with dish soap in hot water.

Cleaning Instructions (Continued)

- Do NOT use solvents or abrasives.
- Rinse with hot, clean water.
- To clean the Electric Pump Main Unit, wipe it with a damp cloth, using only water or water mixed with a mild soap.

Tubing Cleaning Instructions:

- Wipe the outside of the Tubing with a clean, damp cloth. Avoid getting water or other fluids inside the Tubing.
- If liquids are found inside Tubing, rinse the Tubing with clean hot water. Hang to AIR DRY COMPLETELY BEFORE USING. Do NOT insert any object into Tubing.

Note: *Before each use visually inspect the individual components for cracks, chips, tears, discoloration or deterioration; as they may affect the efficiency of the Rumble Tuff® Electric Breast Pump. In the event that damage to the device is observed, please discontinue use until the parts have been replaced.*

Breastmilk Storage Guidelines

Check with your healthcare provider or an International Board Certified Lactation Consultant for specific storage instructions. This information is based on current research and applies to mothers who have healthy, full-term babies who are storing their milk for home use.

Important Notes:

- **STORAGE CONTAINERS:** The best option for storing human milk is glass or hard-sided plastic containers with well-fitting tops. When you are freezing your breastmilk, do not fill the container more than 3/4 full to allow for expansion from freezing. Store milk in the main refrigerator or freezer compartment, away from the door, to avoid changes in temperature that may compromise the milk.
- **RE-WARMING THE FROZEN MILK:** Frozen milk can be thawed in the refrigerator overnight or under cool running water. Gradually increase the temperature of the water to heat the milk to feeding temperature. Be sure to use thawed milk within 24 hours. Discard any leftovers.
- **RE-WARMING THE REFRIGERATED MILK:** Refrigerated milk can be warmed by running under warm water for several minutes. ALWAYS test the temperature of the milk on the inside of your wrist before feeding it to your baby. The milk should feel warm, not hot.
- **WHAT TO DO:** Storing milk in 2 to 4 oz amounts may reduce waste. Write the date you expressed your milk on the container (include your baby's name on the label if your baby is in a daycare setting). Use the oldest milk first. Throw away previously frozen milk that is not used within 24 hours.
- **WHAT NOT TO DO:** Do NOT heat milk directly on the stove. Do NOT thaw in hot or boiling water. Do NOT microwave milk, as microwaving alters the composition of the milk and may create hot spots that may burn your baby. Do NOT refreeze thawed milk.
- **THAT'S WEIRD:** Human milk may separate into a milk layer and a cream layer when it is stored. This is normal. Gently swirl to redistribute the cream before giving to baby. Sometimes thawed milk may smell or taste soapy. This is due to a breakdown in the milk fats. The milk is safe and most babies will still drink it.
- **IS IT DANGEROUS:** Expressed milk can be kept in a common refrigerator at the workplace or day care center. The US Center for Disease Control (CDC) and the US Occupational Safety and Health Administration (OSHA) agree that human milk is not among the body fluids that require special handling or storage in a separate container.

Breastmilk Storage Guidelines (Continued)

Place	Temperature	How long	Things to know
Refrigerator	39°F or colder	<ul style="list-style-type: none"> Up to 72 hours is best. Up to 5-8 days is okay for very clean expressed milk. 	Store milk in the back of the main body of the refrigerator.
Freezer	0°F or colder	<ul style="list-style-type: none"> Up to 6 months is best. Up to 12 months is okay. 	Store milk toward the back of the freezer where temperature is most constant. Milk stored at 0°F or colder is safe for longer durations, but the quality of the milk might not be as high.
Smaller cooler with a blue-ice pack.	59°F	24 hours	Keep ice packs in contact with milk containers at all times; limit opening cooler.
Countertop, table	Room temp (60°F - 85°F)	<ul style="list-style-type: none"> Up to 3-4 hours is best. Up to 6-8 hours is okay for very clean expressed milk. 	Containers should be covered and kept as cool as possible; covering the container with a clean cool towel may keep milk cooler. Throw out any leftover milk within 1 to 2 hours after the baby is finished feeding.

Do not refreeze thawed milk. Discard any remaining thawed milk if not used within 24 hours.

Tips

Pumping Tips:

The most important and successful part of pumping is the let-down reflex. The let-down reflex is when your milk begins to flow freely.

- Try relaxing with a warm, wet compress on your breasts.
- Drink plenty of fluids.
- Avoid outside distractions.
- Visualize yourself happy, well rested and successful at pumping.
- Sometimes a light circular massaging of the breasts prior to and/or during a pumping session will help your let-down reflex start/continue.
- Go to a non-stressful environment.
- Use a comfortable chair that gives lower back support with or without a footstool. After all, you're the one who has to be comfortable.

Depending on your health and milk production condition, you may be able to achieve more than one let-down in one pumping session. Most mothers have more than one let-down while breast feeding and don't even notice the change. To achieve multiple let-downs; which mimics the natural reflex found in breast feeding, you can switch back and forth between the Swift Expression and Natural Expression just by pressing the Power/Mode Button.

Remember that successful pumping is a learned skill. Don't get discouraged, and keep working at it to obtain the best results. Breastmilk is produced on a supply and demand basis; expressing your breastmilk regularly stimulates and maintains your milk supply.

If it feels uncomfortable or painful on your breasts while using the Pump, try decreasing the suction level by pressing the Vacuum Adjusting buttons (-). You may also try pumping with/without the Silicone Cushion. If these tips do not help, turn the product off and remove the pump from your breasts for a few minutes. If you are still experiencing pain or discomfort, consult your healthcare professional or an international board certified lactation consultant.

Tips (Continued)

Battery Storage:

- When storing the battery for a short period of time, store within the following temperature range: -4°F to 140°F (-20°C to 60°C).
- When storing the battery for an extended period of time, store within the following temperature range: 32°F to 77°F (0°C to 25°C).

Breastfeeding Tips:

- The Rumble Tuff® Electric Breast Pump helps you maintain a supply of milk and can help you prepare for planned events such as "Date Night" or unexpected events such as power outages, medical treatment, or relief from tender nipples.

Troubleshooting

Interruption of Electrical Service	<p>If you experience an interruption of electrical service, due to a power surge or other factors, during normal operation of your Rumble Tuff[®] Electric Breast Pump, try the following:</p> <ul style="list-style-type: none"> ● Unplug the AC/DC Adapter cord from the electrical outlet and wait for five seconds. ● Plug the AC/DC Adapter cord back into the electrical outlet. ● If problem persists, try a different electrical outlet.
Radio Frequency Energy	<p>This equipment generates, uses, and can radiate radio frequency energy. If not installed properly and used in accordance with the instructions; it may cause harmful interference with other devices in the vicinity. Try the following:</p> <ul style="list-style-type: none"> ● Reorient or relocate the device receiving interference. ● Increase the distance between the equipment. ● Connect the equipment into an outlet on a different circuit from that which the other device(s) are connected. ● Consult the manufacturer or a field service technician of the device receiving interference for help. ● There is no guarantee that interference will not occur.
Decreased Suction	<p>If it feels like the suction has decreased, try the following:</p> <ul style="list-style-type: none"> ● Make sure the Breast Pump is assembled correctly (see Assembly Instructions pages 8-13). <ul style="list-style-type: none"> ○ Inspect the One-Way Valve for any tears or pinholes. This component is critical to achieving proper suction. If damaged, replace the One-Way Valve with one of the spare valves that is provided with your Rumble Tuff[®] Electric Breast Pump. ● Detach the Tubing and Tubing Adapter from the back of your Electric Breast Pump and place a finger over the Tubing Adapter Inlet. <ul style="list-style-type: none"> ○ If you do feel suction, the Electric Pump Base is operating properly. ○ If you do NOT feel suction or feel that the suction has decreased to an unexpected strength; contact Rumble Tuff Customer Service for warranty service (see Page 32).

Troubleshooting (Continued)

Pumping Doesn't Start	<ul style="list-style-type: none"> ● Make sure the adapter is properly plugged. ● If you are using the battery, make sure the battery is fully charged. ● When battery runs low, Both Blue & White Lights Flash and then turn off automatically.
Breastmilk Buildup in the Breastshield Body	If breastmilk doesn't appear to be draining fast enough into the Collection Bottle, reduce the suction setting on your Electric Breast Pump to allow the milk to drain before gradually increasing.
Where to get Replacement Parts	<p>If you are missing a part or a part was damaged, please contact Rumble Tuff® customer service by calling toll-free (855) 228-8388.</p> <p>You can purchase additional accessories for your Rumble Tuff® Electric Breast Pump by the following methods:</p> <ul style="list-style-type: none"> ● Email: sales@rumbletuff.com ● Fax: (801) 796-2688 ● Toll-Free: (855) 228-8388

Additional Accessories:

Parts	Part Number	Parts	Part Number
● Breastshield Body with Adaptor Ring.....	PA91500-0101	● AC/DC 12V Power Adapter, UL	PA902U1-0101
● Silicone Cushion (2pcs)	PA90300-0201	● Bottle Adapter, Wide	PA900P0-0000
● One-Way Valve(4pcs)	PA90800-0420	● Collection Bottle Combo Kit (2 pcs).....	PA91900-0201
● Bottle, Wide Mouth with Feeding Nipple (2 pcs).....	PA90400-0201		
● Storage Bottle only, Wide Mouth with Cap (6 pcs).....	PA90401-0601		
● Electric Pump Diaphragm(2pcs)	PA91200-0220		
● Electric Pump Diaphragm Kit	PA91300-0101		
● Manual Pump Diaphragm Kit	PA90900-0101		
● Vacuum Tube with adapter, for Double Pump	PA90600-0120		
● Vacuum Tube with adapter, for Single Pump	PA90700-0120		

Technical Data

Power Requirements:

- AC/DC Adapter (included) North American Plug configuration:
Brand: Acute Ideas
Model: GDP12AV-0901300, 100V~240V AC, 50/60 Hz, 0.45A, wall plug type AC/DC Adapter with DC 9V/1.33A output.
- Battery:
Li-Ion Rechargeable battery; 7.4V 850mAh
Brand: Acute Ideas
Model: 18-0740850-RT00, compatible with Nikon EN-EL9
- The power adaptor plug is also served to disconnect the device, not to position the equipment to make it difficult to operate the disconnection device.

Technical Specifications:

VACUUM:

Reflex/Swift Expression: -115~-200mmHg (-15.3~-26.6kpa)




Natural Expression: -115~-250mmHg (-15.3~-33.3kpa)

Ambient Temperature range: +10°C to 40°C (50°F to 104°F)

Relative Humidity range: 20% to 80%

Transportation / Storage Environment : -10 to 50°C / 20 to 90% relative humidity

Key to Box Symbol:

-  Type BF equipment (Breastshield)
-  Consult instructions for use
-  Disposal of Electrical & Electronic Equipment (WEEE)

Made in China

Warranty

Limited Consumer Product Warranty (United States)

This Rumble Tuff® Electric Breast Pump is warranted to the original consumer "the purchaser" to be free from defects in material and workmanship which are not commercially acceptable for the period of one year from the date of purchase on the pump mechanism and a period of 120 days on all other components. Warranty coverage terminates if you sell or otherwise transfer this product to another person. This warranty gives you specific legal rights and you may also have other rights, which vary by location.

RUMBLE TUFF® MAKES NO EXPRESS WARRANTY OF ANY KIND REGARDING THIS PRODUCT OTHER THAN THOSE WARRANTIES SET FORTH HEREIN. ANY IMPLIED WARRANTY, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR ANY IMPLIED WARRANTY OF FITNESS FOR A PARTICULAR PURPOSE TO THE EXTENT PERMITTED BY LAW, SHALL BE LIMITED IN DURATION TO A PERIOD OF 120 DAYS FROM THE DATE OF PURCHASE BY THE ORIGINAL PURCHASER.

In the event that this product is found by Rumble Tuff® to not meet the above limited warranty, as purchaser's sole and exclusive remedy, Rumble Tuff® will repair or at the option of Rumble Tuff®, replace this product without charge for such replacement parts or labor. The purchaser shall bear all expenses related to returning this product to Rumble Tuff®. This warranty does not apply to any part of the product that has been subject to misuse, abuse, or alteration. Only Rumble Tuff® brand replacement parts are recommended for use with Rumble Tuff® breastfeeding products. Improper or incorrectly performed maintenance or repair, or use of non-Rumble Tuff® brand replacement parts voids this warranty.

TO THE EXTENT PERMITTED BY LAW, RUMBLE TUFF® SHALL NOT BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES INCLUDING, BUT NOT LIMITED TO, REPLACEMENT COSTS RESULTING FROM THE BREACH OF ANY WRITTEN OR IMPLIED WARRANTY.

Warranty (Continued)

If you wish to make a claim under this warranty, you will need the following:

- The entire original Rumble Tuff® product
- The original receipt showing date of purchase
- Return Authorization Number (RA#)
- Detailed description of the problem

You, the purchaser, are responsible for all shipping charges.

To obtain an RA# and shipping address: Please contact Rumble Tuff® Customer Service Department by;

Toll-free (855) 228-8388, or

E-mail to customerservice@rumbletuff.com, or

In writing to:

Customer Service Department

Rumble Tuff

865 N 1430 W

Orem UT 84057

Items returned without a RA# will not be processed and will be returned to sender (postage due).



Rumble Tuff®

Rumble Tuff® is owned and operated by
A Cute Baby, Inc.

865 N 1430 W Orem UT 84057-6402

Tel: (801) 609-8168 Fax: (801) 796-2688

www.rumbletuff.com www.acutebaby.com